

## AN ELECTRON OPTICAL ACHROMAT

Stephen W. Carmichael,<sup>1</sup> Mayo Clinic

Spherical and chromatic aberrations have been the bane of optical lenses ever since they were first ground from a piece of glass. As light travels through a convex (converging) lens, the rays at the center of the optical axis are refracted (bent) less than the peripheral rays, so that the central rays are focused behind the peripheral rays. This is the essence of spherical aberration. Light of differing wavelengths (colors) interact differently with the lens so that longer wavelengths (red) are focused behind shorter wavelengths (blue). This is chromatic aberration. In the early days of light microscopy, these two inherent flaws seriously limited the quality of images. Over two centuries ago, it was first appreciated that these aberrations could be corrected by mating a convex lens with an appropriate concave (diverging) lens of glass having different refractive properties. What this accomplished was to introduce aberrations that were equal and opposite, effectively canceling the original aberrations. This combination of convex and concave lenses is known as the achromat, and is present in every optical microscope, except the cheapest toys. Of course, you can appreciate that this review is a gross simplification, as achromatic lenses are more complicated (therefore expensive) than I have described. But the general principles hold.

Electromagnetic and electrostatic lenses behave as a convex lens, converging a beam of electrons toward its focal point. As you know, the resolution is limited by spherical and chromatic aberrations. Spherical aberration follows the same principle as an optical lens, and chromatic aberration varies only by substituting variations in the energies of individual electrons for wavelength. Manufacturers of electron optical equipment have minimized these aberrations by limiting the beam size and narrowing the energy spread in the electron beam, but the lenses are still of the uncorrected type. Obviously, what is needed for further improvement is a divergent lens, but apparently no such thing exists. A solution for this problem was recently demonstrated by Gertrude Rempfer, Denis Desloge, and Walter Skoczylas of Portland State University and Hayes Griffith of the University of Oregon.<sup>1</sup> Their ingenious solution was to insert an electron mirror into the optical path, introducing aberrations that are equal and

opposite.

What is an electron mirror? In this case it's a converging hyperbolic electrostatic mirror. The key element is a conical electrode that creates a retarding electric field and turns the electrons back (reflects) at a curved equipotential surface in front of the electrode. The electrons in the peripheral part of the beam are less strongly focused than more central electrons. This is the opposite of what happens in a lens, thereby correcting spherical aberration. An electron of higher energy penetrates more into the field of the mirror than an electron of lower energy, and therefore is focused more strongly. Again, this is the opposite of what happens in a lens, and this corrects for chromatic aberration.

Introducing a mirror into the optical path meant that a method for separating the incident and returning beams has to be provided. Rempfer *et al.* devised a Y-shaped optical bench with the electron source on one limb of the Y, a separating magnet where the limbs converge, the electron mirror at the base of the Y, and the final image on the other limb. They performed tests with the optical bench by putting a fine copper grid near the final image and casting a shadow pattern. In a remarkable series of illustrations with uncorrected lenses, they showed "pin cushion" and "barrel" distortions of shadows of the grid made with both visible light and electrons. These esthetically-pleasing images were chosen for the cover of the journal, but the astonishing thing about the images was how similar the distortions were, regardless of whether they were formed by light or electrons. Additional experimental demonstrations showed how the electron-generated images could be corrected with the mirror.

Rempfer *et al.* further explain how improvements in resolution resulting from aberration correction will be more readily demonstrated in electron probes and emission microscopes than in transmission electron microscopes because the aberrations are larger in these instruments than in transmission microscopes. However, I would anticipate that this electron optical achromat will prove to be a breakthrough in improving the resolution of all electron optical instruments. I, for one, can't wait to see the improved results! ■

1 The author gratefully acknowledges Dr. Gertrude Rempfer and Dr. Hayes Griffith for reviewing this article

2 Rempfer, G.F., D.M. Desloge, W.P. Skoczylas, and O.H. Griffith, Simultaneous correction of spherical and chromatic aberrations with an electron mirror. An electron optical achromat. *Microscopy and Microanalysis* 3:14-27, 1997.

### Front Page Image

## SEM Photograph of a Black Fly (Simuliidae)

Imaged on an ETEC SEM at 5 kV using the SEM Wideband Multi-Detector Color Synthesizer (designed, built and patented by David Scharf). Then acquired digitally at 2,048 X 1,538 pixels directly into a Macintosh Power PC as a TIFF file, using Digital Micrograph software and Digiscan hardware. Then output to a CELCO film recorder, using Ektachrome 100+ film, to produce a 4 X 5 transparency.

Photograph is Copyright by David Scharf, 1994

All Rights Reserved

As one of a series, this image is available at 350X as a 19" X 23" high quality, lithographic print.

For information on the series, visit the Microscopy Today WWW site at <http://www.microscopy-today.com>

Or visit Microscopy Today at Booth #1017, Microscopy & Microanalysis '97

## MICROSCOPY TODAY

A monthly newsletter dedicated to the unique interests in microscopy and mailed to some 8,000 microscopists worldwide - at no cost in the U.S. and with a modest international charge.

PO Box 620122, Middleton, WI 53562 - Tel.: (608)836-1970 - Fax: (608)836-1969 - eMail: [MicroToday@aol.com](mailto:MicroToday@aol.com)

<http://www.microscopy-today.com>


Copyright in 1997 by Today Enterprises. All rights reserved.

Don Grimes, Editor




# Topcon Introduces Opti-SEM.™

## A Whole New Way Of Looking At An Old Problem.

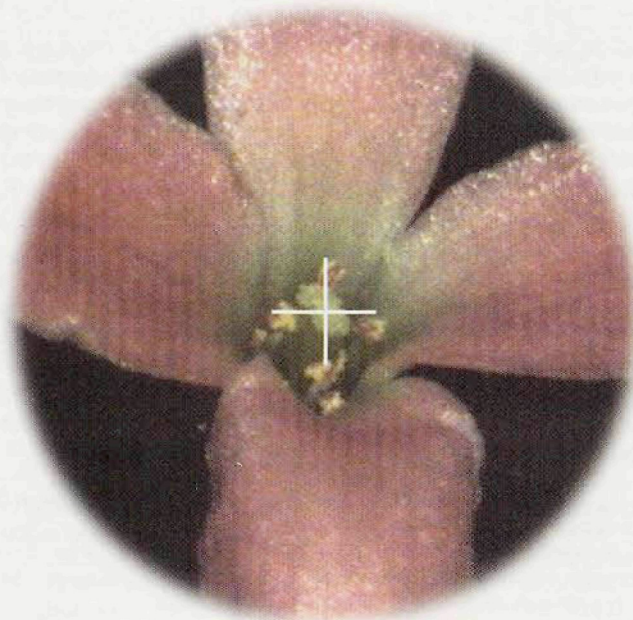

Optical Image

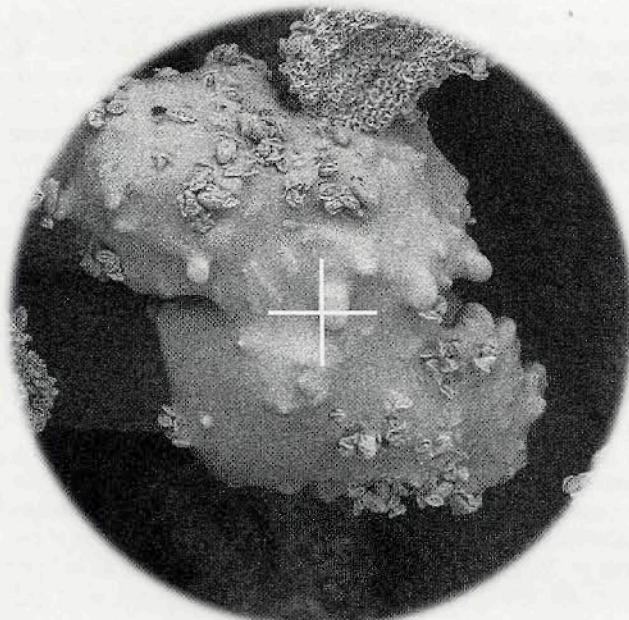

SEM Image

The Built-In Targeting Microscope Lets You Actually Find What You Want To See.  
Think What That Means For Productivity.

With magnification, sometimes more *is* less. Locating a SEM target on a large, complex surface can waste a lot of time. Every time.


#### An Old Astronomer's Trick.

Astronomers use a small, broad field telescope to aim the main telescope. How easy. How productive. So Opti-SEM borrowed the idea. Its built-in 5x to 40x optical microscope makes targeting just as easy for you.

In fact, everything is easy with Opti-SEM. The Windows GUI interface is so familiar, the operator needs little training. And image averaging assures a clean, crisp image every time. No matter whether you output to digital, video or Polaroid.

#### Topcon's #1 Target—Affordable Flexibility.

Opti-SEM's stand-alone, industry-standard computer helps keep you at the state-of-the-art for less. Hardware and software upgrades are local and quick. Now when technology

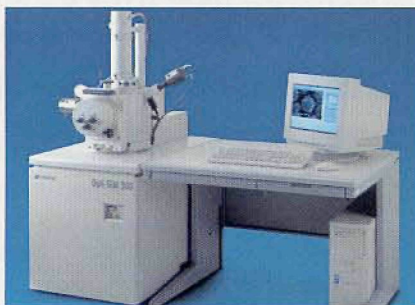

advances, you can afford to move with it.

Better yet, there's an Opti-SEM 300 to meet your particular needs. Including a low-vacuum version. So set your sights on the best value in SEMs today. Topcon. Call 1-800-538-6850 for complete information.



Our Focus Is Clear...To Provide The Best Value  
In SEMs Today.



**LAST CALL!**

Yes, Microscopy & Microanalysis '97, the joint meeting between The Microscopy Society of America, the Microbeam Analysis Society and the Histochemical Society, is about to take place - August 10-14 in Cleveland, Ohio.

OK, so Cleveland might not be the place you would select for your (next) honeymoon. Thanks, however, in major part to the efforts of the local arrangement committee, an excellent program is scheduled - including scientific papers, exhibition of the latest microscopy related products, poster and micrograph competitions and, of course, social events.

Should you be making a late decision to attend, you can find a listing of hotels at the MSA web site: <http://www.msa.microscopy.com>

★ The National Center for Electron Microscopy announces its fourth Annual Summer School on Computer Interactive High-Resolution Transmission Electron Microscopy, including image acquisition, image processing and image simulation, to be held at the National Center for Electron Microscopy during the week of August 25-29, 1997.

The aim of the school is to train participants in the techniques of computer-assisted high-resolution electron microscope image acquisition and image interpretation, including remote-control microscopy. Participants will learn general principles and apply them to specific cases. Participants will be taught the use of computers to obtain images on NCEM microscopes, followed by training in the use of application programs for image interpretation by image processing and image simulation. Participants wanting to apply school techniques to their own projects will be encouraged to extend their visit to NCEM into the next week - note that this requires a proposal be submitted with

advance notice sufficient for project approval.

For information, see: <http://ncem.lbl.gov/NCEM/workshops.html>

★ The June '97 Microscopy Listserv Archives are not on line. You may find them on the MSA WWW site under: <http://www.msa.microscopy.com>

★ **Looking for a new job?**

The Histology Home Page now has a section for listing job openings. Go to: <http://www.vgernet.net/histology/>  
Then, near the bottom of the page, you will find a section for listing jobs. Fill out the form, and your position will be listed for one month - free. And for those looking for a job, look there, also. There is even a section for those looking for temporary help.

★ **A Web Site of Interest!**


The following web site includes a periodic table - plus, where you may click on an element and get all sorts of useful information, like solubilities in different acids: <http://the-tech.mit.edu/Chemicool/>

★ Steven Slap of Energy Beam Sciences, Inc. and Peggy Wenk of the William Beaumont Hospital's School of Histotechnology, have developed a new www site, "The Histotech's Home Page", in conjunction with their "Introduction to the Internet for Histotechs" workshop. The new site is located at: <http://www.vgernet.net/histology/index.html> and has links to many resources relating to histology, light microscopy and immunohistochemistry.

# MRS-3

*a ISO-9000 standard for Microscopy*

Our second generation, traceable, magnification reference standard for all types (SEM, Optical, STM, AFM, etc.) of microscopy. Usable from 10X to 50,000X with pitch patterns from 2µm (±0.1µm) to 8mm (±0.25µm). Pattern height traceable to 0.1mm ± .001µm. New patterns for particle size counting (circles, squares and rectangles) and chemical imaging. Several hundred MRS's have been sold to date.

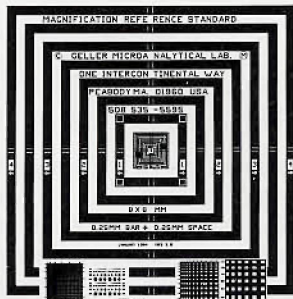

Send for our *free* RESOURCE GUIDE which discusses magnification measurement, error determination and calibration procedures.



426e BOSTON STREET (RT. 1) • TOPSFIELD, MA 01983  
508/887-7000 • FAX: 508/887-6671 • [jg@gellermicro.com](mailto:jg@gellermicro.com)  
<http://www.gellermicro.com>

Certified to the ISO-9001 standard by Lloyd's Registry Quality Assurance  
**"Doing it Right - the First Time"**



**Analytical Services:**

- Electron Probe ● SEM ● Metallography
- Surface Analysis (ESCA & Auger) ● Microhardness

**Products:**

- Elemental Reference Standards (some traceable)
- Magnification Standards (traceable)
- Vacuum Desiccators
- Computer Control System Updates for Electron Probes
- Digital Imaging for Analog SEMs
- Energy-dispersive electronics updates: pulse processor, pulse height analyzers and software (qualitative and quantitative)



426e BOSTON STREET (RT. 1) • TOPSFIELD, MA 01983  
508/887-7000 • FAX: 508/887-6671 • [jg@gellermicro.com](mailto:jg@gellermicro.com)  
<http://www.gellermicro.com>