

MicroscopyEducation


Empowering Pre-College Teachers Through Investigations of Micro- and Nano-Worlds

V. M. Serio, Jr.¹, K. D. Moulton¹*, J. L. Jamison¹, G. Fletcher¹, L. Flower², and S. M. Duboise¹

¹Department of Applied Medical Sciences, University of Southern Maine, 96 Falmouth Street, 178 Science Bldg., Portland, ME 04103

²Farmington River Regional School District, 555 North Maine Road, Otis, MA 01253

* kmoulton@usm.maine.edu

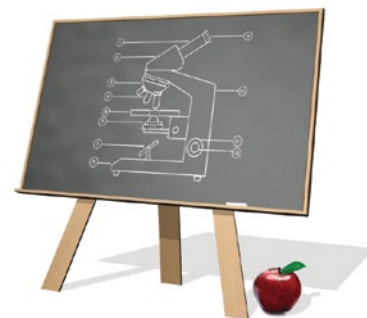

Introduction

The microscopic world provides teachers with excellent opportunities to develop hands-on, inquiry-based curriculum, yet many K-12 educators lack the knowledge, equipment, and skills necessary to engage their students in these types of explorations. Through a Science Education Partnership Award (SEPA) from the National Center for Research Resources (NCRR), researchers at the University of Southern Maine (USM) are offering professional development activities designed to build skills in microbiology, immunology, nanotechnology, and microscopy—subjects that are notably underrepresented in the K-12 curriculum and in the undergraduate education of many teachers [1, 2]. Because these teachers lack formal

microbiology training, the subject matter is often viewed as too dense, abstract, or inaccessible to be presented to pre-college students, despite the immensely important and fundamental roles microbes and viruses play in evolution, human health, and the ecology of all biomes on earth [3-6]. Microscopy is the essential gateway for making the invisible micro- and nano-scale worlds visible, accessible, and less abstract for students at all levels.

Micro- and Nano-Scale Explorations for K-12 Educators

Teacher professional development with the pervasive goal of making microbiology, virology, and micro- and nano-scale imaging accessible for K-12 classrooms is being pursued through two-week summer institutes and weekend academic year programs. In July 2008 the first “Micro- and Nano-space Explorations of Health and Disease” Summer Institute was held at USM. The 14 teachers of grades 3-8 who enrolled were targeted because these grade levels are critical for the development of scientific interests among the students, but teachers at these levels have few classroom microscopy resources and typically have insufficient training in their use. Emphasis during the summer institute was on curriculum development, outreach partnerships, and provision of a light microscope with digital image capture capability in each participant’s classroom. To reveal the nano-scale during the professional development, the project incorporated extensive use of USM’s transmission electron microscope (TEM). Sessions with the teachers involved specimen preparation for TEM (methods described below) and acquisition of images using the scope. Teachers viewed the stomata on lettuce leaf specimens with light microscopy as shown in Figure 1A and chloroplast ultrastructure in a tomato leaf with TEM as shown in Figure 1B. Viruses, particularly bacteriophages, as shown in Figure 1C, were also important specimens for the TEM work. We recently described elsewhere our partnership with high school educators to support classroom use of bacteriophage projects for instruction about viruses and molecular biology [7]. By arrangement with the FEI Company, a benchtop scanning electron microscope (SEM) was used by teachers to acquire images of diatomaceous earth, shown in Figure 1D. Initial teacher responses to both summer and academic year

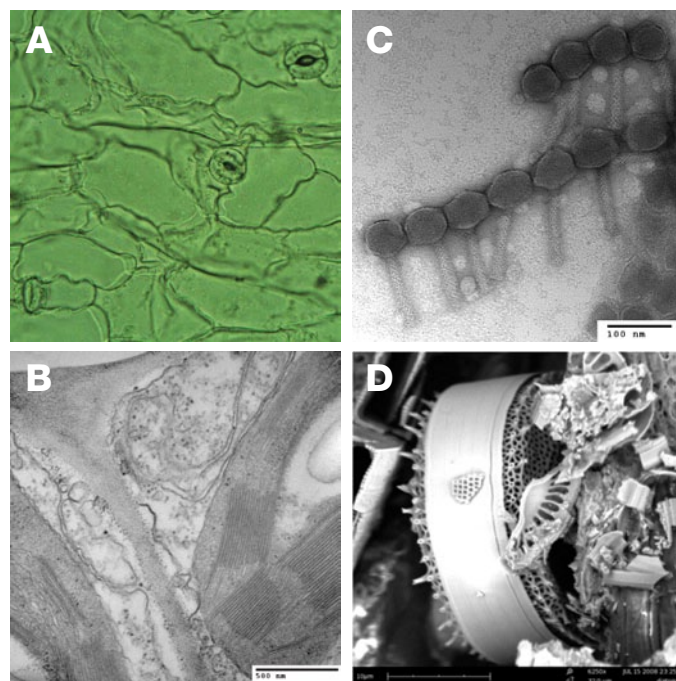


Figure 1: Light and electron micrographs acquired during the USM teacher professional development summer and academic year sessions. (A) Lettuce leaf by light microscopy showing stomata in the open and closed position upon the addition of a drop of saltwater (Image width = 300 μ m). (B) TEM thin section of tomato leaf showing cell wall and chloroplast thylakoid layers. (C) TEM image of negatively stained bacteriophages (that infect *Escherichia coli*, strain C) isolated by teacher participants during the Fall 2008 session. (D) SEM of diatoms. Image acquired using the FEI Phenom during the 2008 teacher summer institute.

Expand your Knowledge of Microscopy with MSA Membership!

Whether your primary focus is in light, electron, or scanning probe microscopy, or the biological or the physical sciences, MSA takes your knowledge to the next level!

Members Receive:

- A personal subscription to MSA's official journal, *Microscopy and Microanalysis*, and MSA's popular bi-monthly magazine, *Microscopy Today*.
- Peer Networking through the Society's Focused Interest Groups and Local Affiliated Societies.
- Discounts on books, journals and other educational materials.
- MSA Awards Programs, Scholarships, Speaker Opportunities, and more!



Join MSA Today!

Visit: www.microscopy.org or call 1-800-538-3672

offerings have been highly positive. A future phase of the project will include the dissemination of project curriculum products and the integration of academic year programs for K-12 teachers into USM's academic offerings.

Materials and Methods For Transmission Electron Microscopy

Plant tissues were fixed for two hours in 3% glutaraldehyde in 0.1 M cacodylate buffer at pH 7.4. During the fixation, nitrogen gas was bubbled into the plant tissue samples. Tissue samples were rinsed with three changes of cold 0.1 M cacodylate—0.1 M sucrose buffer, 15 minutes per change, and post-fixed in cold 2% osmium tetroxide for 2 hours on ice. The specimens were subsequently washed for 15 minutes with 0.1 M cacodylate buffer. Plant tissue samples were then dehydrated through a series of cold ethanol: 30%, 50%, 70%, 80%, 90%, 100%. The 100% ethanol step was followed with one 30-minute change of cold (1:1) solution of 100% ethanol—propylene oxide. The dehydration series was completed with three changes of 100% propylene oxide at 4°C, 30 minutes per change. Plant tissues were infiltrated with a 2:1 mix of Spurr's resin—propylene oxide for 2 hours followed by two changes of 100% Spurr's resin for two hours each. Samples were polymerized at 65°C for at least 12 hours. Sections 60-80 nm thick were cut with an ultramicrotome and poststained with 2% uranyl acetate and lead citrate. Thin sections were observed using a TECNAI G2 Spirit BioTWIN TEM operating at 100 kV. Teacher participants observed the procedures but did not handle toxic reagents.

Outcomes and Future Goals

Each teacher cohort attending the summer sessions acquires more advanced skills in light microscopy and better understanding of the various types of microscopy. Consequently, they are able to make better use of existing microscopy resources in their schools and to take advantage of the full range of capabilities of their new microscope with built-in digital camera. Educators exit the workshop thinking more broadly about the applications of microscopy. Units in mathematics, life and physical sciences, and literacy created over the summer indicate that teachers see microscopes as tools to engage learners across the curriculum. After some initial classroom testing, these lessons will be made available to teachers worldwide on the supporting website. This program has extended and enhanced the educational outreach of our virology laboratory and TEM facility. Teacher workshop participants have been engaging researchers at our institution to advise on curriculum development, to present talks to their students, and to provide student tours of university laboratory and microscopy facilities. We look forward to augmenting our workshops to further inspire teachers' creative use of microscopy in their curriculum and to encourage regional pre-college educators to participate in additional formal university science courses that are increasingly being offered at times that permit in-service teachers to attend. Response to the project work indicates ample interest of educators in developing their scientific knowledge and skills and in discovering how to better incorporate microscopy into their curriculum.

Acknowledgments

This work was supported by Award Number R25RR024280 from the National Center For Research Resources, a component of the National Institutes of Health. The content is solely the responsibility of the authors and does not necessarily represent the official views of NCRR or NIH. TEM resources have combined support from National Science Foundation (grant CNS-0521262) and the University of Southern Maine. Assistance from FEI representatives Vanessa Mazuzan, Diane K. Stewart, and Cliff Mathisen with the FEI Phenom SEM is gratefully acknowledged. [MT](#)

References

- [1] M Schaechter et al., *Microbiology in the 21st Century: Where we are and where are we going*. American Academy of Microbiology, Washington, DC, 2004.
- [2] M Bloom, *Frontiers in Microbiology*, BSCS, Colorado Springs, CO, 2006.
- [3] JA Fuhrman, *Nature* 399 (1999) 541.
- [4] KE Wommack and RR Colwell, *Microbiol Mol Biol Rev* 64 (2000) 69.
- [5] CA Suttle, *Nature* 437 (2005) 356.
- [6] M Schaechter et al., *Microbe*, ASM Press, Washington, DC, 2006.
- [7] SM Duboise et al., *The Science Teacher* 76 (2009) 32.

**Why have over 2,000 scientists
in 35 countries selected
Minus K[®] vibration isolators?**

**Our Negative Stiffness systems
deliver 10x to 100x better
performance than air systems and
even better than active systems.**



Without Minus K[®]



With Minus K[®]

**The best performance *and* the lowest price.
That's hard to beat!**

VIBRATION ISOLATION BY:

minus k[®] TECHNOLOGY

460 S. Hindry Ave., Unit C, Inglewood, CA 90301
Tel: 310-348-9656 Fax: 310-348-9638
sales@minusk.com • www.minusk.com



You asked. We listened.

Ultrahigh Resolution

Analytical Performance

Ultralow kV

Uncoated stretched polypropylene

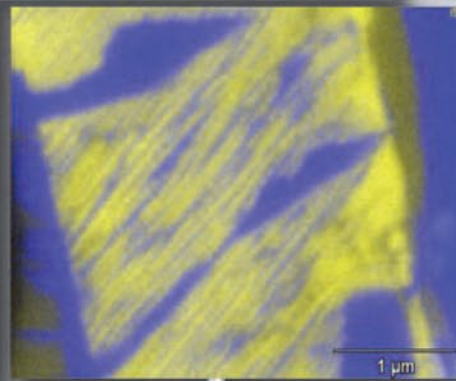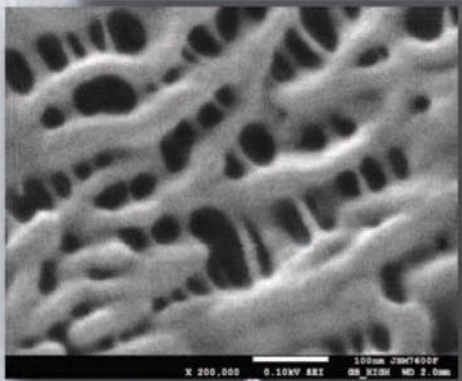


EDS/WDS

EDS map of TiO₂ in FeO_x with < 100nm spatial resolution 30,000x

LABe

1-3 nm twinning in mineral BSE image


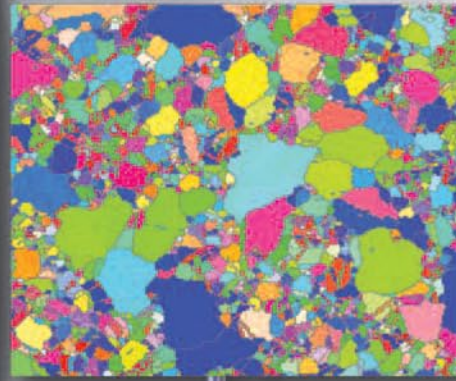


EBSD

Orientation map of Ni alloy

STEM

CNT with 1-3 nm Pt nanoparticles

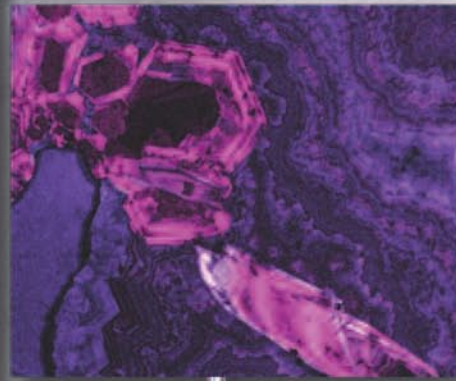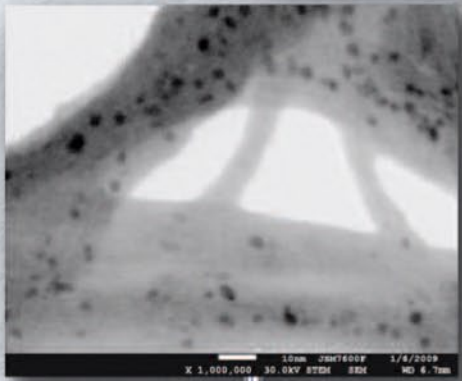

CL

Quartz and zinc oxide

Seeing is believing!

Contact salesinfo@jeol.com to arrange a demo.
For more information or to view a brochure visit www.jeolusa.com/7600F.

All magnifications as printed at 4 x 5 in.



JSM-7600F

Another
Extreme
Imaging
Solution

JEOL

Stability • Performance • Productivity
www.jeolusa.com • salesinfo@jeol.com
978-535-5900

Thermal FEG SEM