


Microscopy TODAY

2021 Buyers' Guide



Achieve cutting edge results with DiATOME...

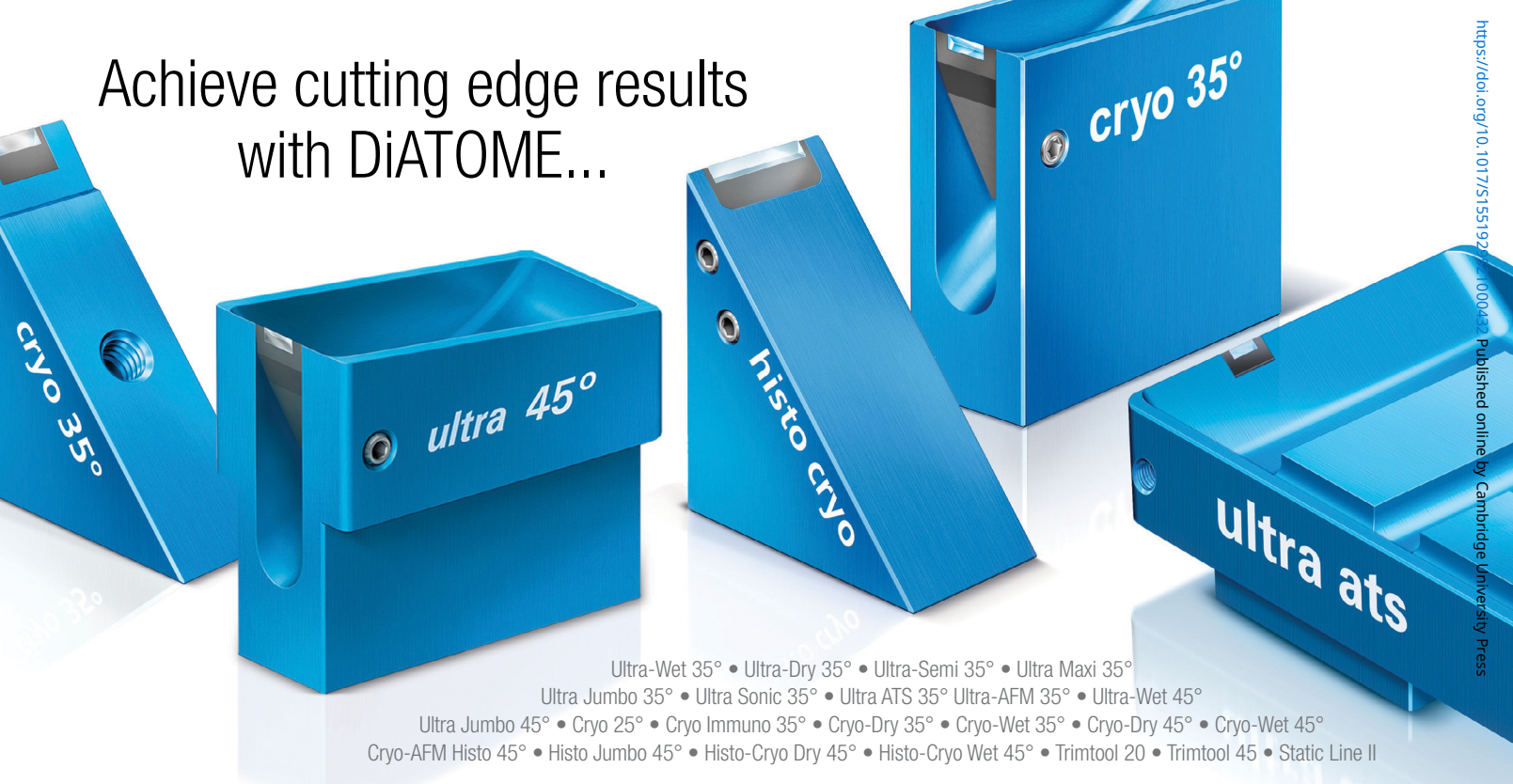

Innovation, using DiATOME Diamond Knives...

Micro-Optical Sectioning Tomography to Obtain a High-Resolution Atlas of the Mouse Brain

Existing imaging tools have limitations for brainwide mapping of neural circuits at a mesoscale level. In collaboration with DiATOME, researchers developed a Micro-Optical Sectioning Tomography (MOST) system utilizing a DiATOME Diamond Knife that can provide micron tomography of a centimeter-sized whole mouse brain.


Slicing was performed by moving the specimen to generate ribbons, and each ribbon was simultaneously imaged. The illuminating beam passed through a beam splitter, mirror and objective to irradiate the ribbon. The imaging beam collected by the objective and passed through the mirror, beam splitter and tube lens was then recorded by a line-scan CCD.

A 3D structural dataset of a Golgi-stained whole mouse brain at the neurite level was obtained. The morphology and spatial locations of neurons and traces of neurites were clearly distinguished. Researchers found that neighboring Purkinje cells were sticking to each other.

Acknowledgement

Micro-Optical Sectioning Tomography to Obtain a High-Resolution Atlas of the Mouse Brain Anan Li, Hui Gong, Bin Zhang, Qingdi Wang, Cheng Yan, Jingpeng Wu, Qian Liu, Shaoqun Zeng, Qingming Luo

Britton Chance Center for Biomedical Photonics, Wuhan National Laboratory for Optoelectronics—Huazhong University of Science and Technology, Wuhan 430074, P. R. China.

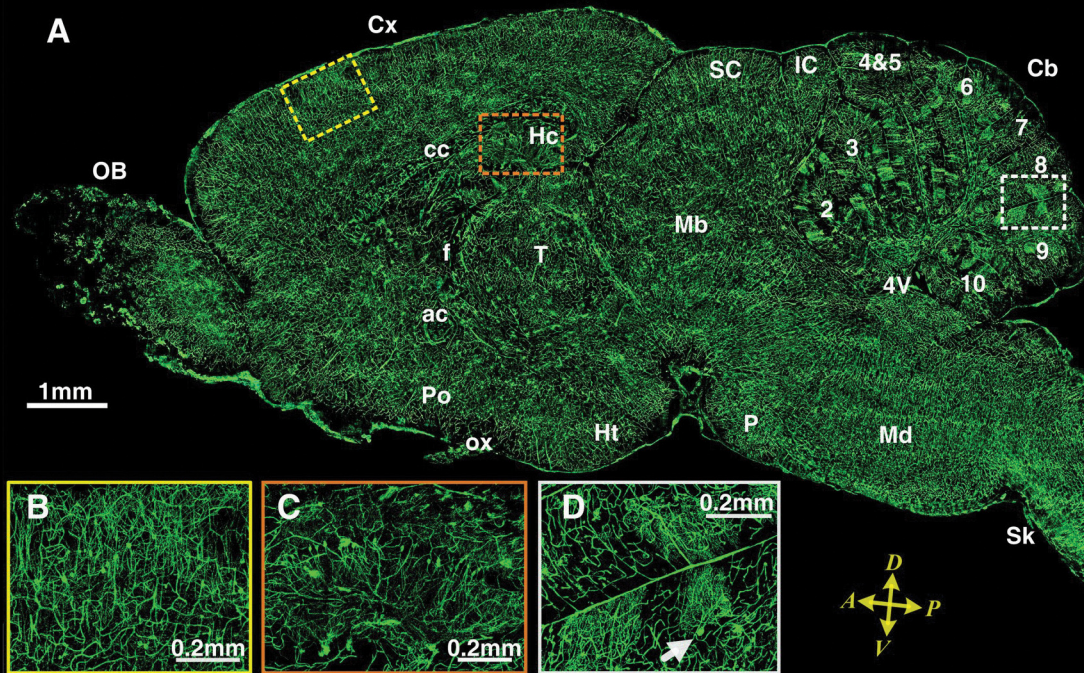

DiATOME U.S.

for more information, please visit our website at...

www.emsdiasum.com

P.O. Box 550 • 1560 Industry Rd.
Hatfield, PA 19440

Tel: (215) 412-8390

Fax: (215) 412-8450

email: info@emsdiasum.com

or stacie@ems-secure.com

Welcome to the Microscopy TODAY Buyers' Guide.

The guide is published in print and online at www.cambridge.org/mto. For such a guide to be truly useful, it must be comprehensive. Thus, the editors have assembled an index identifying hundreds of products listed under more than 300 categories.

HOW TO USE THIS GUIDE:

Print version: (1) Look up the product category in the Products and Services Index.
(2) View the product descriptions in the Microscopy Product Vendors section.
(3) View more complete product descriptions at the company website.
(4) View the scope of certain companies in the Company Profile section.

Online version: (1) Go to www.cambridge.org/mto.
(2) Click on "2021 Buyers' Guide."
(3) View an exact replica of the Buyers' Guide and use the search box in the upper right to find a particular kind of product.
(4) Use the Company Index section to find a particular company (page numbers here refer to Company Profiles and Major Listings).

To companies not listed in this guide:

To list your company and promote your products in our next edition, please contact Kelly Miller at kmiller@mrvica.com or 856-768-9360.

This publication is a supplement to *Microscopy Today*. It is published by Cambridge University Press for the Microscopy Society of America, a not-for-profit society dedicated to the promotion of microscopy and microanalysis.





Hybrid-pixel electron detectors for Electron Microscopy

- The highest count rate combined with high frame rate
- Unprecedented dynamic range
- Noise-free detection

Find out more at www.dectris.com/EM